


FlexNeedle™

Sample Like Never Before

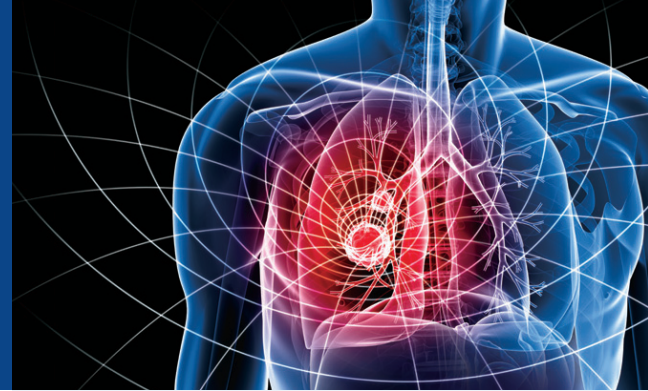

FlexNeedle's innovative flexible shaft and proprietary coring tip enable cutting and capturing of core specimens for histology and cytology.

Unique Advantages

- Innovative flexible needle with coring tip
- Prevents damage to the bronchoscope during insertion and withdrawal by locking the needle tip inside a flexible coiled-steel protector
- Needle internal lumen is 0.95 mm enabling passage of a stylet (included) or other devices



Easy-to-Use Handle

- Ergonomic, single-handed needle actuation
- Adjustable depth stops provide penetration control from 5 to 20 mm
- Standard luer fitting is compatible with wide range of aspiration syringes and stopcocks (not included)

FlexNeedle™ – Disposable, single-patient-use (box of 5)

| Catalog Number | Needle Gauge | Penetration Depth (mm) | Maximum Catheter OD (mm) | Catheter Length (mm) |
|----------------|--------------|------------------------|--------------------------|----------------------|
| 10005 | 18 | 0–20 | 1.9 | 1400 |



Access targets in hard-to-reach locations with the proprietary flexible needle shaft

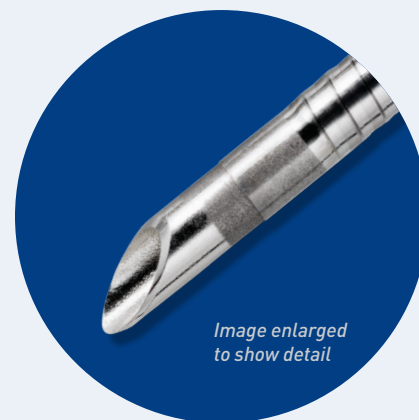


Image enlarged to show detail

5 mm coring tip is laser-welded to the metal coil shaft

Order Information:
877.428.1600
sales@broncus.com

Copyright © 2010 Broncus Technologies, Inc. Indications, Contraindications, Warnings and Instructions for Use can be found in the product labeling supplied with each system.
CAUTION: Federal (USA) law restricts this device to sale by or on the order of a physician. Rx only.

MK-271 rev 8/10